

Made in Korea

Innovation of Respiratory Care

OmniOx


Multi-Functional Non-Invasive Respiratory Care System

MEK

Innovation of Respiratory Care


OmniOx

OmniOx, A Fully Featured Ventilator

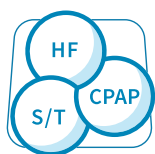

A Fully Featured Ventilator

Provides extensive oxygen therapy from HF to NIV

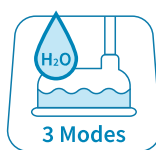

Adaptive Control Humidifier

Support optimal humidity by adaptive control algorithm



Visible Graphic Wave

PTIF, ROX Index, Flow, FiO2
PIP, PEEP, Vtest, RR, SpO2, PR



CMS & Wifi & Remote Control

Enable immediate clinical response without expose to infection

OmniOx, New Standard of HFNC

Minimize medical staff effort

PTIF(Peak Tidal Inspiratory Flow) Viewer

- Display the actual PTIF in real-time
- Ensure optimal flow rates

ROX Index

Assist in deciding the continuation of therapy

Continuous patient therapy

Bi-Flow Modes

Reduce respiratory resistance and increase efficiency by adjusting inhalation and exhalation flow rates differently

TSF

(Pulse oximeter with target SpO2 feedback control by FiO2)

- Maximize time within the target SpO2 range
- Reduce modulation of SpO2 and workload the caregiver



Minimizing Medical Staff's Effort



Wave Monitoring System

Visualization of parameters for setting up oxygen therapy

Optimal flow rate setting

- Real-time visualization of patient's breathing status
- Ability to provide flow rate exceeding the patient's maximum inspiratory flow


Appropriate assessment of patient's therapy

- Possible to monitor cannula fitting status and changes

Prevention of ROX Index intubation failure

- Real-time monitoring of optimal indicator for evaluating HFNC results

● Bi-Flow

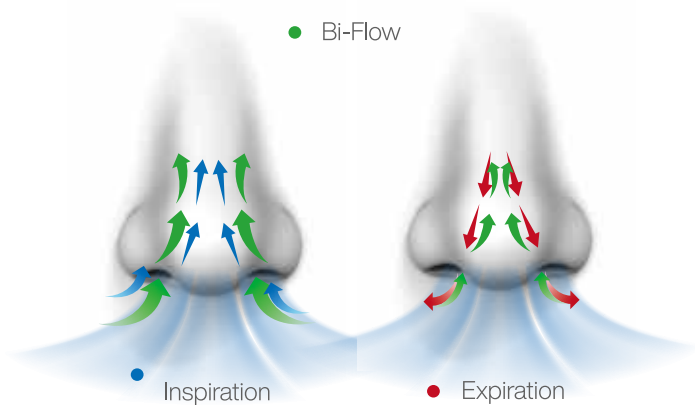

Bi-Flow Mode

Improvement in respiratory efficiency and enhancement of device adaptability


During Inhalation

- Providing a high-flow oxygen
- Effectively removing dead space of nasal cavity
- Preserving the advantages of the high flow nasal cannula

During Exhalation

- Providing a relatively low flow of oxygen
- Reducing expiratory breathing resistance
- Addressing the limitation of the high flow nasal cannula

Continuous Patient Therapy

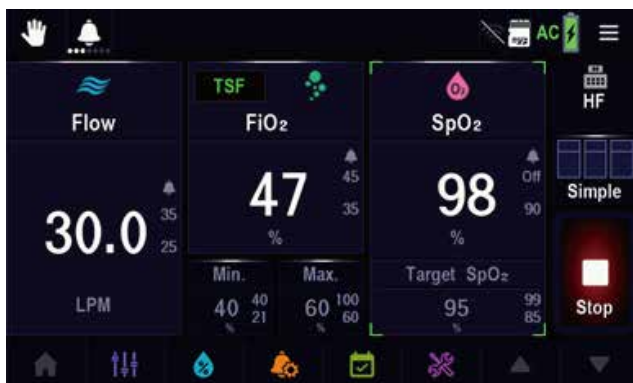

TSF Pulse oximeter with target Spo2 feedback control by FiO₂

- Maximize time within the target SpO₂ range
- Reduce modulation of SpO₂ and workload the caregiver

With Battery

- Provide heated and humidified gas consistently even during intra-hospital transfers

Specificaion

HF (High Flow)

Flow (Inspiratory flow rate) 1 ~ 60 LPM

CPAP (Continuous Positive Airway Pressure)





CPAP 4 ~ 20 cmH2O ±(1.7 + 4 % of the set value) cmH2O
 Pressure Assist + (PA(+)) OFF, 1 - 3 cmH2O
 Pressure Assist - (PA(-)) OFF, -3 - -1 cmH2O
 Trigger Level 3 - 20 lpm
 Rise time Fast(0.2 s), Medium(0.3 s), Slow(0.4 s)
 Auto Start ON, OFF

bi-level Positive Airway Pressure (Spontaneous, Spontaneous / Timed, Timed)

IPAP 4 - 40 cmH2O ±(1.7 + 4 % of the set value) cmH2O
 EPAP 4 - 20 cmH2O ±(1.7 + 4 % of the set value) cmH2O
 Pressure Assist - (PA(-)) OFF, -3 - -1 cmH2O
 Respiration Rate 2 - 60 bpm
 Inspiratory Time(Ti) 0.3 - 3.0 seconds
 Trigger Level 3 - 20 lpm
 Trigger Type S/T, T, S
 Rise time Fast(0.2 s), Medium(0.3 s), Slow(0.4 s)
 Auto Start ON, OFF

- Pulse oximeter can be provided with the two options: MEKICS, Masimo

Product Line

	OmniOx 700	OmniOx 750A	OmniOx 750	OmniOx 751
	Hospital	Hospital	Hospital	Home
Color				
Features				
Battery	-	✓	✓	✓
Modes	HF	HF	HF, NIV(CPAP,bi-level)	HF, NIV(CPAP,bi-level)
Functions				
TSF	✓	✓	✓	-
Bi-Flow	-	✓	✓	✓
PTIF* Viewer	✓	✓	✓	✓
ROX Index	✓	✓	✓	✓
Cannula Fitting	✓	✓	✓	✓