


# Introducing S-scan the Next Step in Office MRI

Based on extensive customer feedback and years of engineering, **Esaote** has designed the **S-scan with eXP Technology**, an optimized MRI scanner for any practice with a substantial musculoskeletal work-load. The **S-scan** covers all musculoskeletal anatomy, from foot to shoulders including the most important spine segments including L and C-spine.

## S-scan from Esaote, quality in MSK imaging

Since 1993, **Esaote** has been developing and producing MRI systems specifically for musculoskeletal imaging. Esaote MRI systems are designed and optimized specifically to meet the demanding needs of MSK professionals.

S-scan integrates our 20+ years of MRI experience into a system that provides you outstanding image quality and an easy to use and install system.

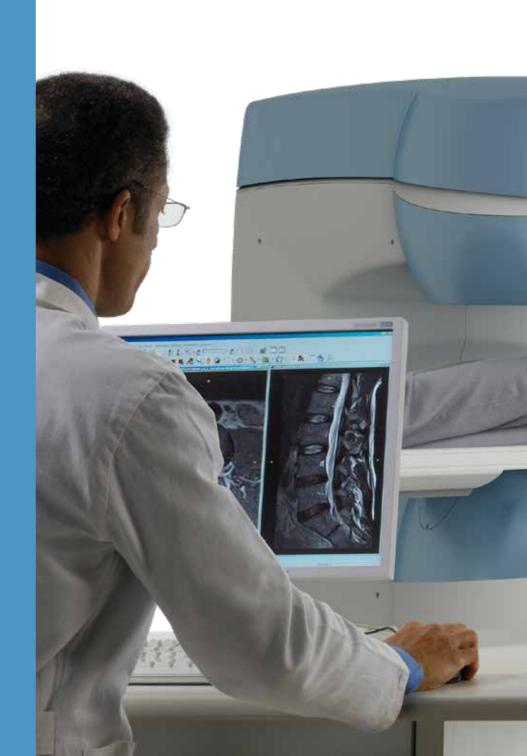







## **Esaote market leader** in MSK MRI systems

Musculoskeletal MRI is the second most requested MRI study behind neuro imaging. Esaote's focus is MSK MRI and we have become the world leader in that area.

Esaote systems are used daily in many clinics and hospitals all over the world for routine clinical use as well as for research. Esaote systems are unique as they are extremely easy to use and combine excellent image quality with optimum patient comfort.

## Efficiency with S-scan musculoskeletal MRI

An efficient MRI system will provide you: ease of use, simple patient workflow, cost effectiveness and outstanding image quality. S-scan has been developed to provide all of those features and more. The open design for fast and easy patient positioning, gantry display with real-time, DICOM worklist,

network integration, no cryogens and everything in a small foot print. It all adds up to make S-scan an efficient and cost effective MRI regardless of whether you perform one or twenty examinations per day.


#### **Esaote Dedicated MRI**

Esaote's ever-growing family of innovative dedicated systems have been on the market for 20 years and over 2400 systems have been sold worldwide.



### **@XP** Technology

Technology provides faster overall system response thanks to a number of sophisticated MRI software techniques and a high-power workstation.

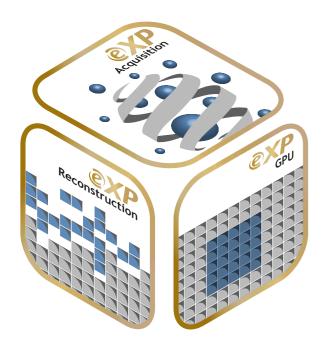

















#### **EXP** Performance enhancements

By combining sophisticated acquisition and reconstruction methods such as SpeedUp and TR reduction (Patent pending), eXP will substantially improve the MRI examination by reducing scan times and improving image quality.







## S-scan patient positioning

#### S-scan: easy patient access

Easy patient access is one of the features that distinguishes the S-scan as an ideal system for more advanced applications such as MRI arthrography. The wide opening and compact magnet provide easy access to the patient.

## S-scan: open gantry and rotating table

S-scan features a very wide, asymmetric rotating table. This wide table not only makes patient positioning very easy, it also gives maximum stability and comfort to the patient for an easy setup and better images.

#### S-scan: real-time positioning

Fast error free positioning with the real-time feature that includes fast sequences and a display on the system gantry. Start scout scan directly from the gantry and check patient positioning without running back and forth to the console. The simplest and fastest way to check correct patient positioning.

#### S-scan coils

The S-scan comes with a complete set of RF coils, from foot to spine.



Shoulder Coil DPA



Knee Coil DPA



Wrist/Hand Coil DPA



Ankle/Foot Coil DPA



Optional
Cervical Spine
Coil DPA



Optional 4 Channel Lumbar Spine



Optional Lumbar Spine Coil DPA



Cervical Spine
Coil
Complete with
flexion and extension
options



**Optional Flexible**Coil



Optional Bllateral TMJ Coil

### S-scan: connected to MSK needs

#### **Connectivity**

The S-scan includes complete network, archiving and documentation features to work either stand alone or as part of an integrated environment.

As a stand-alone scanner: S-scan comes with an integrated DVD archive and retrieve software package, printer output facility and a patient-CD package.

As a networked system: S-scan is DICOM compliant, and offers smart solutions for connectivity to your PACS and teleradiology.

#### Windows® interface

The S-scan is easy to learn as it uses Windows® functionalities.

The interface and protocols are custom designed for musculoskeletal imaging which speeds up and simplifies

the examination procedure considerably.

#### **Remote Service**

All Esaote systems come with remote service capability. Many system parameters and system components can be checked via the remote service program specially developed for Esaote MRI systems. The remote

service capability allows
Esaote to react to any
potential issue and to
get you up and running
faster!



#### S-scan economics

S-scan is one of the most economical MRI systems available, featuring easy installation, the patented Esaote RF pavilion, low energy consumption, no cryogens, and remote service capability. These features allow S-scan to be suitable for sites with a limited patient work-load.


#### S-scan easy to install and maintain

S-scan is a "one room" MRI system which means that the complete system, magnet, electronics and console can be installed in a single room of only 20 m² (238 ft²).


- S-scan is based on an optimized permanent magnet so no cryogens and no complicated cooling systems are required.
- Dedicated shielding available. The Esaote dedicated shielding is a pavilion style independent shielding that can be installed without any renovation to the space.
- Fast and high quality service. Thanks to the built-in remote service capability, technical assistance is fast and efficient.

#### S-scan, compact configuration

The S-scan also offers a compact configuration (optional) for smaller installation sites. A specially designed compact bed is provided which can be accommodated in an RF pavilion of just 9 m<sup>2</sup> (100 ft<sup>2</sup>).

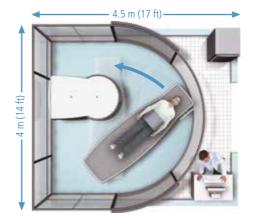

**Standard Configuration** 

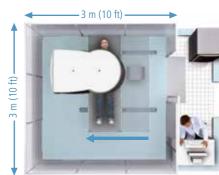

**Compact Configuration** 









ESAOTE S.p.A.- sole-shareholder company
Via Enrico Melen, 77 16152 Genova, ITALY, Tel. +39 010 6547 1, Fax +39 010 6547 275, info@esaote.com

